



etherea^{MX}

#1 LASER PLATFORM
IN THE WORLD'S 2ND
LARGEST AESTHETIC
MARKET.



#1 LASER PLATFORM
IN THE WORLD'S 2ND
LARGEST AESTHETIC
MARKET.

70+

TREATMENT INDICATIONS
AND GROWING

**SUPERIOR
SERVICES**
TO MEET
YOUR NEEDS

**2-YEAR
WARRANTY**
OPERATE WITH
CONFIDENCE AND
PROTECT YOUR
BUSINESS

A VERSATILE
SYSTEM FOR
THE **GLOBAL
MARKET**

**MARKET
LEADER**
OVER A 2,000
SYSTEMS
INSTALLED
WORLDWIDE

**MULTIPLE
HANDPIECES**
WITH MULTIPLE
OPTIONS



#1 LASER PLATFORM IN THE WORLD'S 2ND LARGEST AESTHETIC MARKET.

- A **LASER** platform with **6 technologies** and **70+ tx indications**.
- **Market share leader** in the 2nd largest aesthetic market worldwide
- **Maximum ROI** in an all-in-one system with numerous available handpieces.
- **2-year** warranty
- **Number one** in Brazil with more systems sold than all others combined.




ETHEREA® was born in Brazil and expanded globally around the world in just over 13 years, spreading to Europe, Asia, Oceania and North America. Since 2009, ETHEREA® has continued to grow, incorporating business models that have made the platform a remarkable sales success: from the ETHEREA® Users' Meeting aimed at developing and expanding clinical practice with technology, to after-sales assistance and service, creating a solid foundation for any medical practice.



QUALITY,
PERFORMANCE
AND LONGEVITY



A NEXT GENERATION SYSTEM:
With additional technologies and expanded indications, ETHAREA-MX® is a valuable asset for any practice.

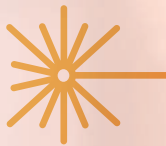

MANUFACTURED TO GLOBAL STANDARDS:
Manufactured to ISO/UL standards and having obtained FDA, CE, CAN/CSA and ANVISA CERTIFICATIONS, ETHAREA-MX® is a platform suited for the global market.



ESTABLISHED TECHNOLOGIES, CUTTING EDGE EXECUTION:
Efficient fractional delivery, reducing residual thermal damage. Square-wave pulse technology that provides safe and improved outcomes. Versatile treatment options from an all-in-one plug and play handpiece system. These are just some examples of the ETHAREA-MX® advantage.

TECHNICAL SPECIFICATIONS

	ETHAREA-MX
TECHNOLOGY	LIGHT and LASER platform
INPUT VOLTAGE	Automatic Dual Voltage (110-220V)
DIMENTIONS (L x H x D)	48 x 88 x 43 cm
WEIGHT	52.5 Kg
PLATFORM SYSTEM:	YES
Nº OF HANDPIECES AVAILABLE:	6

GREATER VERSATILITY WITH UNIQUE ALL-IN-ONE FILTER AND SPOT TECHNOLOGY.



A single platform with
over **70 treatment**
indications**

DEVICE	AESTHETIC INDICATIONS*	MEDICAL INDICATIONS
 IPL-Sq®	• Pigmented lesions • Scars and striae • Vascular lesions • Hair removal	• Acne vulgaris • Rosacea
 ACROMA-QS®	• Tattoo removal • Non-ablative skin resurfacing (LASER toning)	• Melasma • Pigmented lesions (lentigos)
 LongPulse®	• Vascular lesions • Scars and striae • Non-ablative skin resurfacing (thermal peeling) • LASER hair removal	• Rosacea • Onychomycosis
 ProDeep®	• Non-ablative fractional skin resurfacing • Scars and striae • Wrinkles	• Acne vulgaris
 GoSmooth®	• Non-ablative fractional skin resurfacing • Scars and striae • Wrinkles	• Melasma
 DualMode®	• Fractional ablative skin resurfacing • Scars and striae • Wrinkles • Pigmented lesions	• Vaginal atrophy

*This product is certified as a medical device in the European Union under the Medical Device Directive 93/42/EEC by POLISH CENTRE FOR TESTING AND CERTIFICATION CE 1434, No 1434-MDD-422/2021, exclusively for the indication(s) declared as medical. Other nonmedical uses ascribed to this device are not within the scope of CE certification, and users should be aware product performance and/or safety has not been evaluated by POLISH CENTRE FOR TESTING AND CERTIFICATION for the aesthetic indications herein stated.

**The number of total indications include medical and aesthetic indications, as well as the subtypes of lesions.

IPL-Sq®

- **Proprietary SQUARE-WAVE PULSE®** technology with micro processed pulses for enhanced safety and efficacy;
- **All-in-one system** with 6 available cut filters at 400, 515, 540, 580, 640 and 695 nm, considerably increasing the versatility and profitability of the handpiece;
- **More practicality with cut filters** in the plug-and-play model, with automatic recognition;
- **Special VascuTip® spots** in all wavelengths, working with the utmost efficiency in the treatment of localized lesions and less accessible regions;
- **Built-in sapphire contact cooling system**, to enhance patient comfort and procedure safety;
- 5 to 100 ms pulses for **safety in the treatment of darker phototypes** (Fitzpatrick IV-V), and **effective outcomes for vascular lesions** of different sizes and depths;
- **Repetition rates of up to 2 Hz and a broad treatment area**, measuring 4.8 cm², making procedures much quicker;



All-in-one cut-off filters;
400, 515, 540, 580, 640 and 695 nm



TECHNICAL SPECIFICATIONS

	IPL-SQ®
IPL WAVELENGTH	400-1,200 nm
CUT-OFF FILTERS AVAILABLE	400, 515, 540, 580, 640 and 695 nm
SPOT SIZE	40 x 12 mm
VASCUTIP™	Ø 8 mm and 12 x 12 mm
MAXIMUM ENERGY	33 J/cm ²
PULSE WIDTH	5 to 100 ms
COOLING SYSTEM	built-in sapphire contact-cooling
REPETITION RATE	up to 2 Hz

ACROMA-QS®

- **ACROMA-QS® technology is safe for treatment of all phototypes without restriction by time of year**, which allows you to offer your patients treatments even at the height of summer.
- **An excellent investment in Q-switched LASER technology**, adding new possibilities to your clinical practice and opening up a range of new options for cosmetic and therapeutic treatments;
- **ACROMA-QS® for having two wavelengths** - 1064 and 532 - increases the possibilities of treatments being indicated for the removal of tattoos of various colors and all shades of pigmentary lesions.



TECHNICAL SPECIFICATIONS

	ACROMA-QS®
LASER WAVELENGTH	Nd:YAG 1,064 nm
	KTP 532 nm
AVAILABLE SPOT SIZES	1064 nm 3, 5, & 7 mm
	1064 nm 9 mm, 100 mtz/cm² * 532 nm, with 3 & 5 mm
MAXIMUM ENERGY	1,500 mJ
PULSE WIDTH	20 ns
REPETITION RATE	up to 5 Hz



Available spots:
@1064 nm: 3, 5 and 7 mm;
@532 nm: 3 mm;
optional spots:
@1064 nm: fractional 9 mm
and 100 mtz/cm²;
@532 nm: 5 mm.

LongPulse®

- **Due to greater penetration of the wavelength and absorption by hemoglobin**, LongPulse® LASER offers the utmost precision, efficacy and safety for the treatment of vascular lesions in veins and vessels of different sizes and depths;
- **The true LongPulse® technology, with longer pulse times of up to 60 ms**, results in prolonged heating of the treatment target point, increasing the vessel's internal temperature and improving the response;
- **With spot sizes of 2, 3, 6 and 9 mm in diameter and maximum fluency of up to 500 J/cm²**, LongPulse® offers maximum range for the treatment of various medical indications;
- **Epidermal cooling coupler system:** Increased precision, comfort and safety for the patient;
- **With variable pulses of 0.3 to 1 ms, DYNAMICS®** technology offer effectively adding a technological solution to expand the range of clinical indications to include non-ablative skin resurfacing (thermal peeling), rosacea and onychomycosis.



Available spots: 2, 3, 6 and 9 mm.

TECHNICAL SPECIFICATIONS

	LongPulse®
LASER WAVELENGTH	Nd:YAG 1,064 nm
AVAILABLE SPOTS SIZES	2, 3, 6 & 9 mm
MAXIMUM ENERGY	500 J/cm ²
PULSE WIDTH	0.3 to 60 ms
INTEGRATED COUPLING PLUG SYSTEM	Air-cooling sytem adapter available
MODE OF OPERATION	LONGPULSE, DYNAMICS®
REPETITION RATE	up to 10 Hz

ProDeep®

- **ProDeep®: LASER effect with greater depth of penetration when** compared to other non-ablative LASERs available on the market, ensuring treatments with much greater clinical efficacy.
- **Wide variation in treatment intensity with ProDeep®,** combining higher intradermal heat in subdermal tissue, pulse durations in the range of 3-20 ms, micro-fractionated lenses and collimated beams, providing a unique treatment option;
- **ProDeep® has revolutionized the market,** by offering a non-ablative fractional LASER option with excellent results from only the first sessions, while still being a highly safe technology with virtually no downtime. A great option for body procedures;
- **Greater safety, efficacy and usability with a new high quality optical fractionation system,** square spot formats, smaller micro-fractionated beams and greater spacing between spots;
- **The legacy of a technology brought by ProDeep®,** a LASER technology widely used and approved by renowned doctors, with numerous success stories that prove its effectiveness.



TECHNICAL SPECIFICATIONS

	ProDeep®
LASER WAVELENGTH	Nd:YAP 1,340 nm
AVAILABLE SPOT SIZES	6 mm 8 mm, 100 mtz/cm² 10 mm, 400 mtz/cm²
MAXIMUM ENERGY	200 mJ/mtz
PULSE WIDTH	3 to 20 ms
INTEGRATED COUPLING PLUG SYSTEM	Air-cooling sytem adapter available
REPETITION RATE	up to 2,5 Hz



available spots: fractional 8 mm and 100 mtz/cm², fractional 10 mm and 400 mtz/cm², and 6 mm;

GoSmooth®

- GoSmooth® incorporates **the world's gold standard in skin resurfacing technology** and adds a new treatment proposition to ETHEREA-MX;
- **True non-ablative fractional resurfacing** with GoSmooth® acts with precision on the epithelial tissue layer while generating a homogeneous and residual intradermal thermal effect;
- **Prevents, maintains and treats all aging signs**, and is indicated for the treatment of wrinkles, scars, melasma, and striae.
- Especially designed to **ensure therapeutic efficacy and safety in procedures for any region of the body**;
- **Provides visible results from the first session** and offers a safe procedure with great clinical efficacy for all skin types, without disrupting daily activities;
- **High-quality beam micro-fractionation** through an improved precision lens system, ensuring a much safer use of GoSmooth® in practice;
- GoSmooth® places new therapeutic possibilities, at your disposal, a **high-value technology with great versatility and real results**.



TECHNICAL SPECIFICATIONS

	GoSmooth®
LASER WAVELENGTH	Er:GLASS 1,540 nm
AVAILABLE SPOT SIZES	8 mm, 100 mtz/cm ² fractional lens array
	10 mm, 400 mtz/cm ² fractional lens array
MAXIMUM ENERGY	95 mJ/mtz
PULSE WIDTH	10 and 15 ms
INTEGRATED COUPLING PLUG SYSTEM	air-cooling sytem adapter available
REPETITION RATE	up to 1 Hz



Available spots:
fractional 8 mm and
100 mtz/cm², and
fractional 10 mm and
400 mtz/cm²

DualMode®



- **DualMode® features** an exclusive dual pulse system (stacking) with vaporization and coagulation effects in the same shot, bringing the real power of fractional resurfacing to clinical practice;
- **Cold erbium ablative effect superficial and controlled:** High LASER absorption by water at a wavelength of 2940 nm, combined with purely ablative shots (300-500 µs), results in controlled surface treatments with minimal downtime;
- **Compared to CO2 LASER, DualMode®** offers treatments with less risk of complications and adverse effects, in addition to reduced downtime and rapid post-treatment improvement.
- **DualMode® allows a combination of treatments with different levels of aggressiveness**, using fractional lenses at 100 and 400 mtz with different diameters, widely variable pulse durations – 300 µs to 5 ms – and repetition rates of up to 5 Hz;
- **DualMode® InLift®:** an intraoral treatment that is complementary to ablative fractional rejuvenation, practically painless, without downtime, producing excellent results.
- **ATHENA®** for treatment and care of women´s intimate health is a revolutionary method that promotes the improvement of women´s intimate health.

TECHNICAL SPECIFICATIONS

	DualMode®
LASER WAVELENGTH	Er:YAG 2,940 nm
AVAILABLE SPOT SIZE	Ø 8 mm, 400 mtz/cm² Ø 8 mm, 100 mtz/cm² Ø 2.5* and 6 mm InLift® for intraoral LASER treatment
AVAILABLE ACCESSORIES	ATHENA® 90, ATHENA® 360
MAXIMUM ENERGY	8 J
PULSE WIDTH	0.3 to 400 ms
INTEGRATED COUPLING PLUG SYSTEM	Smoke evacuator or air-cooling sytem adapter available
REPETITION RATE	up to 5 Hz

*optional spot size

Available spots:
InLift® spot for intraoral treatment, plus fractional 8 mm and 100 mtz/ cm², fractional 8 mm and 400 mtz/ cm², and 6 mm; optional spots: 2,5 mm.



ATHENA® 90+ and ATHENA® 360 for women's intimate LASER treatment.



VYDENCE® Medical
Manufacturing Plant and R&D
Rua Aldo Germano Klein, 359. CEAT
13573-470 São Carlos, SP. Brazil

VYDENCE® Medical
Training Center
Av. Fagundes Filho, 486 SL. 93
São Judas 04304-000 - São Paulo - SP

VYDENCE® Medical
U.S. Office
601 Brickell Key Dr, STE 901
Miami, FL 33131 - USA